





IN16B Reduction with Mantid

Progress

#	Description	Status
PR #16620	Renamed entry in nexus file [technical]	Merged
PR #16707	ILL: Load Multiple Files [technical] (there is no prefix)	Merged
Issue #16753	IndirectILLReduction Algorithm Refactoring* + Handle multiple files + Treatment of active single detectors in the LoadILLIndirect + Control mode → Unmirror options (Verena) → Adjust IndirectCalibration accordingly (see discussion by Verena)	In progress
Issue #16754	IndirectILLReduction GUI* + Allow multiple file selection + Add a checkbox for summing of files + Add a checkbox for Control mode + Add a dropdown for unmirror option → Add corresponding handlers for these new GUI items	In progress (will merge with above)
	ShiftCircular (needed for unmirroring) - A new generic algorithm to perform circular shift of x-axis of the workspace	TBC
	Generalisation of FindEPP (needed for unmirroring) - Remove unneeded x-axis unit validation, update the category	Green light by the author


Single Detectors Treatment

- There are 8 slots in nexus files
- But only few are enabled
- Adjusted the loader

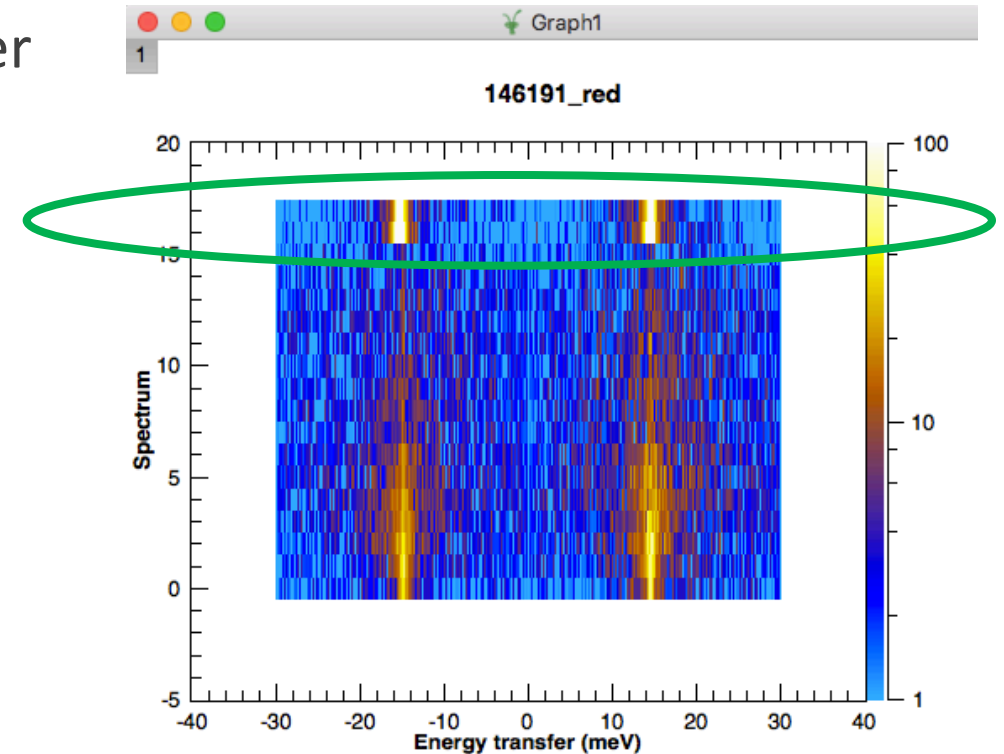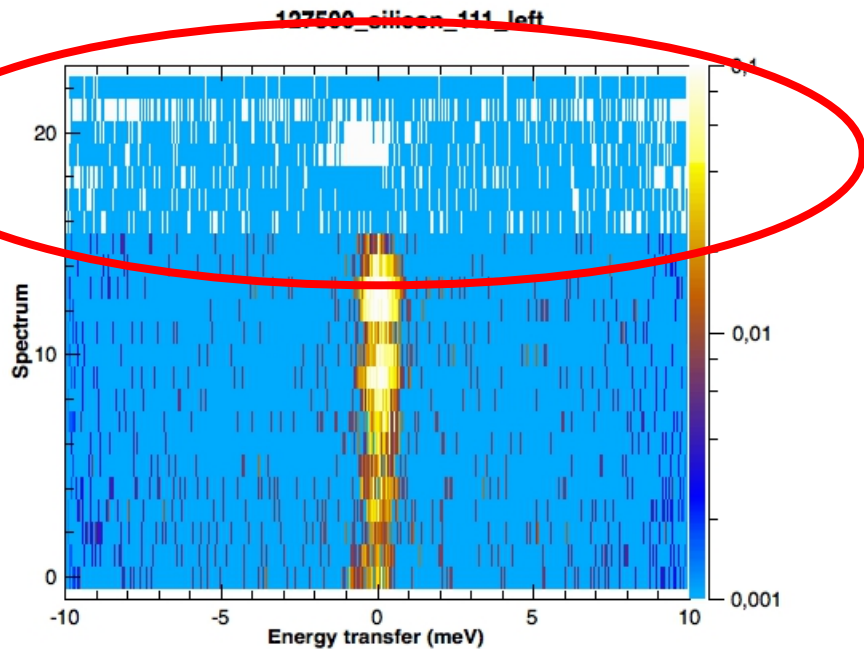

Caveat

- ❑ Loader now does checks from .nxs file which SDs are active, and loads only them
- ❑ But the reduction relies also on Detector Grouping map file (.xml)
- ❑ This file needs to be consistent with the nexus files that one is reducing
- ❑ Not a big problem, just need to be documented
- ❑ What to put in the default .xml file?

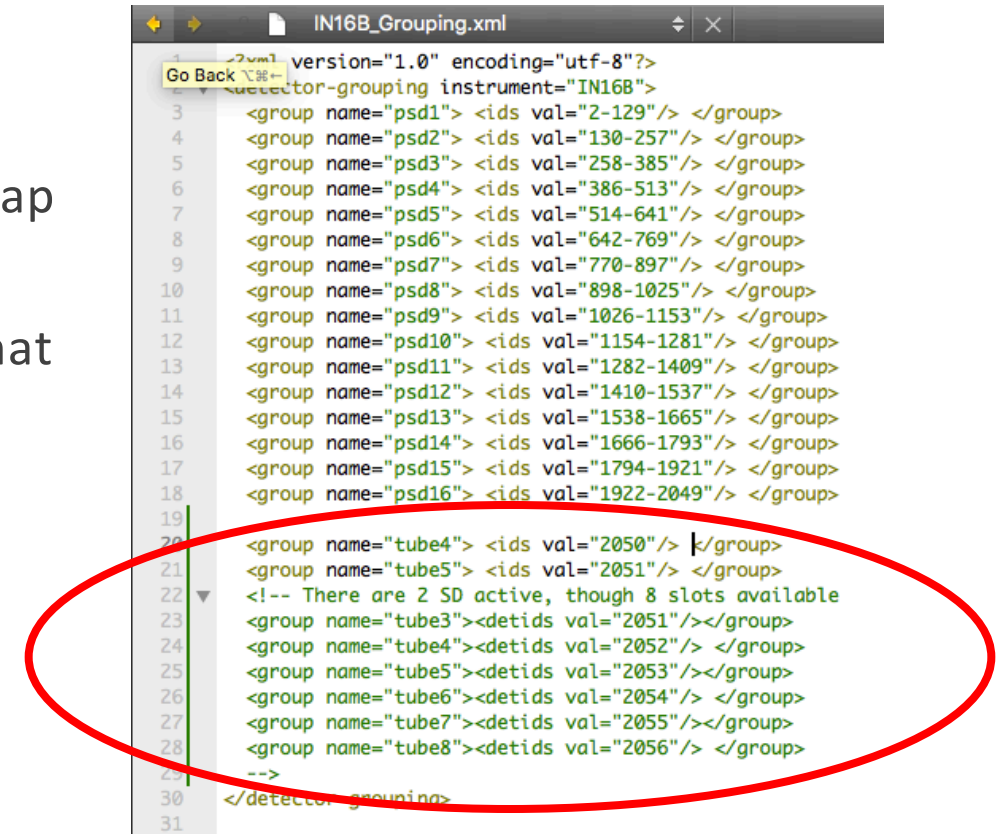

```
1 <?xml version="1.0" encoding="utf-8"?>
2 <detector-grouping instrument="IN16B">
3   <group name="psd1"> <ids val="2-129"/> </group>
4   <group name="psd2"> <ids val="130-257"/> </group>
5   <group name="psd3"> <ids val="258-385"/> </group>
6   <group name="psd4"> <ids val="386-513"/> </group>
7   <group name="psd5"> <ids val="514-641"/> </group>
8   <group name="psd6"> <ids val="642-769"/> </group>
9   <group name="psd7"> <ids val="770-897"/> </group>
10  <group name="psd8"> <ids val="898-1025"/> </group>
11  <group name="psd9"> <ids val="1026-1153"/> </group>
12  <group name="psd10"> <ids val="1154-1281"/> </group>
13  <group name="psd11"> <ids val="1282-1409"/> </group>
14  <group name="psd12"> <ids val="1410-1537"/> </group>
15  <group name="psd13"> <ids val="1538-1665"/> </group>
16  <group name="psd14"> <ids val="1666-1793"/> </group>
17  <group name="psd15"> <ids val="1794-1921"/> </group>
18  <group name="psd16"> <ids val="1922-2049"/> </group>
19
20  <group name="tube4"> <ids val="2050"/> </group>
21  <group name="tube5"> <ids val="2051"/> </group>
22  <!-- There are 2 SD active, though 8 slots available
23  <group name="tube3"><detids val="2051"/></group>
24  <group name="tube4"><detids val="2052"/> </group>
25  <group name="tube5"><detids val="2053"/></group>
26  <group name="tube6"><detids val="2054"/> </group>
27  <group name="tube7"><detids val="2055"/></group>
28  <group name="tube8"><detids val="2056"/> </group>
29  <!-->
30 </detector-grouping>
31
```